



CoolRepair

Scientific UK Ltd

UniFreez™ UND80

Under-Bench SMART-LAB™ -86°C ULT Freezer

Innovative technology to improve the quality for laboratory.

 **DAIHAN**
Scientific
YOUR PARTNER *in* LABORATORY



CONTENTS

UniFreez™ UND80 Features	04
UniFreez™ Freezing System	06
Powered by the Smart-Lab™	08
Green Technology	10
Product Specification	11
About Cool Repair Scientific UK	14
Complementary Services	17





UNIFREEZ™ UND80 UNDER-BENCH SMART-LAB™ -86°C ULT FREEZER



**Innovative
Freezing System**

**SMART
Connectivity**



**Green, Energy
Efficient
Technology**



PRODUCT FEATURES:


- * High quality insulation panel and inner door.
- * Filter free mechanisms, lowering compressor generated heat.
- * Automatic vacuum breaker, enabling easy door re-opening.
- * Self-diagnostic functionality, including built-in warnings and alarms.
- * Automatic temperature data recording.
- * Built-in storage and inventory management.
- * Ergonomic, 7" full touch screen LCD controller.
- * Robust, flexible door handle design.
- * CE Certified.

UNIFREEZ™ FREEZING SYSTEM


INNOVATION AT ITS BEST

DAIHAN Scientific ULT freezers are unique in the way they are built.

Whilst most alternative ULT freezers are designed using a typical cascaded compressor system that employs two compressors working simultaneously, those built by DAIHAN Scientific use a single set compressor instead. Specifically, the 'Simplified Hi-Efficiency Freezing System'.

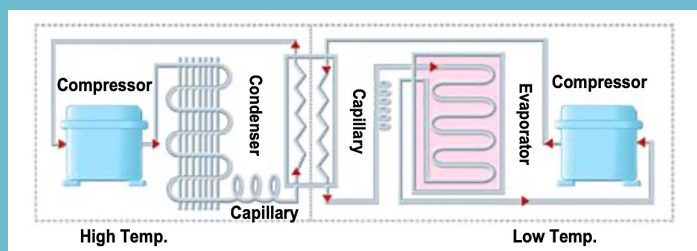

Typical Cascaded Compressor System

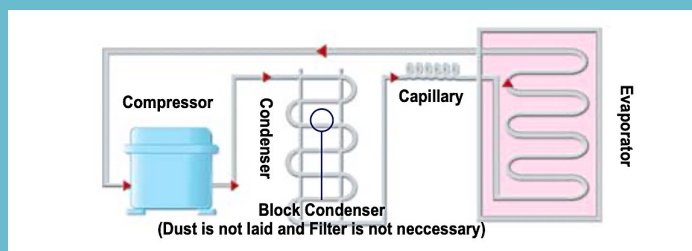

Simplified Hi-Efficiency Freezing System

As illustrated above, the overall freezing mechanism for DAIHAN Scientific ULT freezers is far less complex than that of the typical cascaded compressor system.

In fact, by using this Simplified Hi-Efficiency Freezing System in all ULT freezer models, DAIHAN Scientific overcomes the technical complexity of the typical cascaded compressor systems and simplifies the freezing mechanism all together. Critically, with no negative impact on overall performance.

They also drive real value to end user by:

- reducing energy consumption;
- minimising noise pollution;
- improving equipment sustainability and operational lifetimes;
- increasing efficiency;
- ensuring user safety; and,
- mitigating the risk of operating failures.

DOUBLE THE SECURITY AND TWICE THE RELIABILITY

DAIHAN Scientific offers two unique ULT freezer models, both employing the Simplified Hi-Efficiency Freezing System overleaf. They are:

- **The UniFreez™ ULT Freezer**

The most innovative freezing system in the world. It employs one single Simplified Hi-Efficiency Freezing System to maintain an internal ULT freezer temperature of -86°C .


This is the Freezing System employed by this model of ULT Freezer.

- **The DuoFreez™ ULT Freezer**

An extremely stable freezing system boasting two independent Simplified Hi-Efficiency Freezing Systems which enables an internal ULT freezer temperature of -90°C .



► Comparison Data of Freezing Speed between competing products

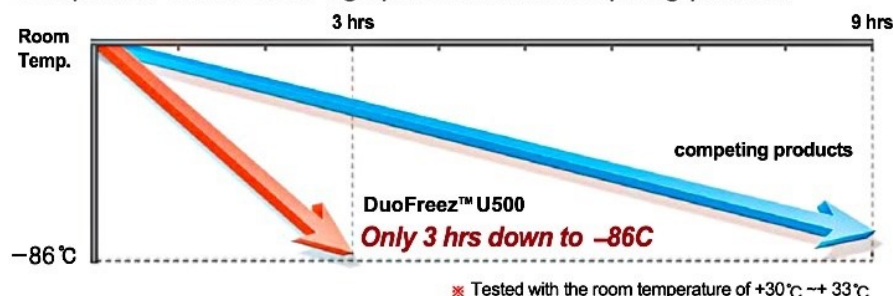

► Time Taken to Reach Down to Ultra Low Temperature

Model	DuoFreez™			UniFreez™				
	U300	U500	U700	U300	U500	U700	WUF-41	WUF-51
Temp.	308Lit.	503Lit.	683Lit.	308Lit.	503Lit.	82Lit.	375Lit.	438Lit.
	Upright Type			Chest Type				
-80°C				240min.	330min.	150min.	120min.	150min.
-86°C	160min.	180min.		270min.	420min.		210min.	270min.
-90°C			240min.			270min.		
-95°C	270min.	300min.						

※ Tested with the room temperature of $30 \sim 33^{\circ}\text{C}$

POWERED BY THE SMART-LAB™

The UniFreez™ U500 -86°C ULT Freezer is continually evolving, like your smart phone in hand and smart TV at home.

All UniFreez™ models come equipped with a 7-inch, full-touch screen LCD controller, integrated Smart-Lab™ System and internet connectivity as standard.

Unique and innovative to the entire DAIHAN Scientific ULT freezer range, this progressive system provides a simple, easy to control, secure and convenient solution to self-diagnostics, real-time monitoring and data recording. Storage and inventory management too.

Using a smart phone, tablet or laptop/computer, users can easily connect their chosen device to the UniFreez™ U500 via the WiRe™ App and/or Web Service.


The Smart-Lab System™ provides real time, immediate access to your equipment - anywhere, anytime.

NOTE: Smart-Lab™ functionality is subject to availability. For further details, speak to a member of our team today.

8 SETTINGS

Set up and adjust various Functions for the unit.


9

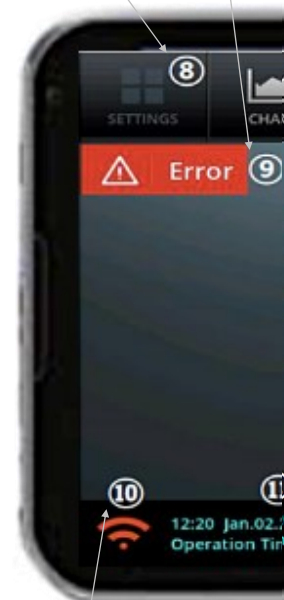

10 WIFI STA

Displays the connectivity

WiRe™ App & Web Service

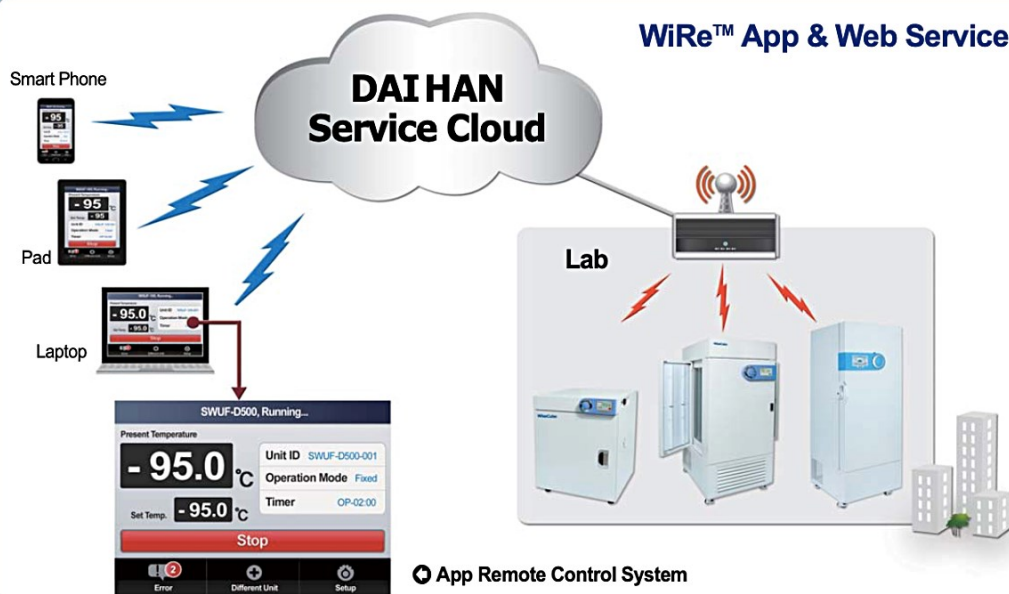

7 CHART

Monitor the operating status via the chart.

4 LCO₂

Prepare for a disaster or emergency by injecting LCO₂.

6 OFF SET

Set the calibration value for the temperature.

3 START/STOP

Adjust the start and stop operation.

5 STORAGE

Enter storage locations for inventory management.

2 CURRENT TEMP

Display the current operating temperature.

ERROR

Notifications in the event of an error.



12 OPERATIONAL INFO.

Shows the state of operation via the LED lamp.

1 SET TEMPERATURE

Check and adjust the temperature setting.

13 ON/OFF ICON

Displays the operating status of the unit.

11 DATE

Displays the current date and time.

TUS

State of internet

GREEN TECHNOLOGY

RoHS COMPLIANT

All UniFreez™ UND80 Smart-Lab™ -86°C ULT Freezers are RoHS Certified.

The RoHS (Restriction of Hazardous Substances) is the directive on the restriction of the use of certain hazardous substances (Pb, Cd, Hg, Hex-Cr, PBB, PBDE) commonly used in electronic and electronic Equipment (2011/65/EU).

39% COST & ENERGY SAVINGS

The Simplified Hi-Efficiency Freezing System, unique to DAIHAN Scientific, consumes less power and drives greater cost savings than other typical, -86°C ULT freezers.

Consumption for a Typical -86°C Freezer*

Power Consumption: 6,570kWh / year
CO₂ Output: 3,212.07 kgCO₂ / year
Electricity Cost: USD 1,051.20

Consumption for a UniFreez™ U500 Freezer*

Power Consumption: 4,015kWh / year
CO₂ Output: 1,962.93 kgCO₂ / year
Electricity Cost: USD 642.4

* Calculated by the electricity-specific emission factors sourced from the IEA (2011a).



PRODUCT SPECIFICATION

FUNCTIONALITY

	UND80 Under-Bench SWUF00082 UNI
High / Low Temperature	Yes
Remote Alarm System	Yes
Alerts for:	
Power Failure	Yes
Sensor Errors	Yes
Temp. Deviation / Error	Yes
Door Ajar	Yes
Automatic Vacuum Breaker	Yes
CE Certified	Yes

ADDITIONAL FEATURES

	UND80 Under-Bench SWUF00082 UNI
Caster Wheels	Yes
External Sensor Hole	Yes
Inner Door	Yes
Temperature Recorder	Yes
Racks & Boxes	Optional

TECHNICAL DATA

	UND80 Under-Bench SWUF00082 UNI
Capacity	82
Cabinet Type	Under-Bench
Freezing System	UniFreez™
Compressor	Hermetic Type
Condenser Type	Block (No Pin, Filter Free)
Refrigerant	Natural (CFC Free)

PERFORMANCE

	UND80 Under-Bench SWUF00082 UNI
Cooling Performance	-86°C
Temperature Range	-65°C to -86°C
Temperature Sensor	PT 100

CONTROL

	UND80 Under-Bench SWUF00082 UNI
Controller	Smart-Lab™
Display	7" Full Touch Screen LCD

DIMENSIONS

	UND80 Under-Bench SWUF00082 UNI
Net Weight (kg)	60
Gross Weight (k)	102
Interior Dimensions (WxDxH)	520 x 485 x 340 mm Material: Stainless-Steel
Exterior Dimensions (WxDxH)	1,065 x 722 x 612 mm Material: Power Coated Steel
Packing Dimensions (WxDxH)	1,350 x 760 x 960 mm

ELECTICAL DATA

	UND80 Under-Bench SWUF00082 UNI
Power Supply (V/Hz)	220-240 / 50

POWER CONSUMPTION

	UND80 Under-Bench SWUF00082 UNI
Watts / Hour	321

Cool Repair Scientific UK

ABOUT US

TRUSTED BY:





The University of
Nottingham

UNITED KINGDOM • CHINA • MALAYSIA

HOSPITAL
INNOVATIONS

ICR The Institute of
Cancer Research



GlaxoSmithKline

LEADING UK REFRIGERATION SPECIALIST

With over 30 years' experience, we are trusted by academic institutions, research facilities, hospitals, clinics, pharmaceutical companies, clinical trial centres and biobanks (to name a few) across the UK.

Renowned for delivering exceptional customer service and a responsive, reliable, trustworthy and effective service, trust us when we say:

"You are in safe hands with Cool Repair Scientific UK as we draw on a wealth of knowledge and experience to deliver exactly what you need, when you need it."

A TRUE ONE-STOP-SHOP

We support the entire lifecycle of laboratory refrigeration equipment, offering a true 'one stop shop' for all your laboratory needs. Not only do we supply an extensive range of refrigeration equipment, from +4°C fridges to low temperature -40°C, -80°C ULT and even -190°C liquid nitrogen freezers, we offer a suite of complementary services too.

To that end, we hold ISO9001:2008 accreditation as well as an Human Tissue Authority (HTA) licence through our sister company, MTS Cryo Stores UK, for the storage and transport of human application and research material, licence number: 22499. For further details, feel free to get in touch.

PURPOSE-BUILT FACILITIES

Based in Nottingham, we are well positioned to arrange the distribution and supply of refrigeration equipment, as well as supporting services, for clients across the UK. Our purpose-built storage facilities boast:

- 24:7 security monitoring and call-our service by specialist, local refrigeration engineers;
- backup freezers and standby units in the event of equipment failure;
- permanent diesel generator in the event of a failed power supply;
- continuous temperature monitoring and alarm systems;
- equipment calibration with UKAS traceable probes; and,
- electronic sample inventory management system.

COMPLEMENTARY SERVICES



SERVICING, MAINTENANCE & CALIBRATION

Proven to extend the lifespan and improve performance of refrigeration equipment. We offer a variety of servicing and maintenance contracts, covering everything from the basics to the most comprehensive of equipment health checks.



BREAKDOWN COVER

A strong measure for any contingency plan. Our breakdown cover acts as a form of insurance, or extended warranty, against future repairs and callout charges. Covering a wide variety of manufacturers and models, we provide cover for equipment up to 15 years old.



EQUIPMENT HIRE & TRANSPORTATION

An effective solution to short and long-term projects, temporary overfills, freezer defrosting and equipment. Using our mini-fleet of specialist transport vehicles, we can easily deliver our hire equipment to locations across the UK. Alternatively, we are able to arrange prompt collection and/or delivery of your very own equipment,.



DISASTER RECOVERY

Critical to safeguarding the contents of any refrigeration equipment, our disaster recovery service is designed to protect you and your material in the unfortunate event of a disaster. Highly flexible, our solutions are tailored to you and your requirements. We can help you protect anything from a single box of samples to the contents of entire fridges/freezers.



Cool Repair Scientific UK

Unit D Marlborough Court
Blenheim Industrial Estate
Bennerley Road, Bulwell
Nottingham NG6 8UY.

+44 (0) 115 975 3743
info@coolrepair.co.uk
www.CoolRepair.co.uk

